



# S9 Escape Auto™

SLEEP THERAPY DEVICE

Part of the S9 Series, S9 Escape Auto provides automatic positive airway pressure therapy and comes packed with technological advances to deliver more ... more performance and more comfort.

**Easier patient acceptance of therapy**  
—engaging design  
Award-winning design combining easy functionality with elegance in the bedroom

**Easy-to-read**  
—colour LCD screen  
Menus and treatment options clearly displayed

**Intuitive**  
—control dial  
Simple and logical menu selection

**No rainout**  
—Climate Control  
Avoiding rainout without compromising humidification




**Convenient follow-up**  
—SD card  
Helps physicians to obtain summary usage data

**Barely there**  
—quiet performance  
78% quieter\* via conducted noise through the tube and 34% quieter\*\* via radiated noise at the device

**Ultimate patient comfort**  
—Enhanced Easy-Breathe™ motor plus EPR  
Provides the smoothest pressure changes and eases the work of breathing to promote compliance


**Optimal therapy delivered automatically**  
—APAP algorithm  
Automatically adjusts pressure to deliver effective and comfortable therapy

**Portable and lightweight**  
—compact design  
Small system with its own discreet travel bag

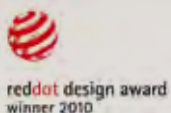

\* Compared to S8™ II devices \*\* Based on ResMed internal studies comparing S9 with devices from a leading competitor

**S9™** Series  
More. Comfort.

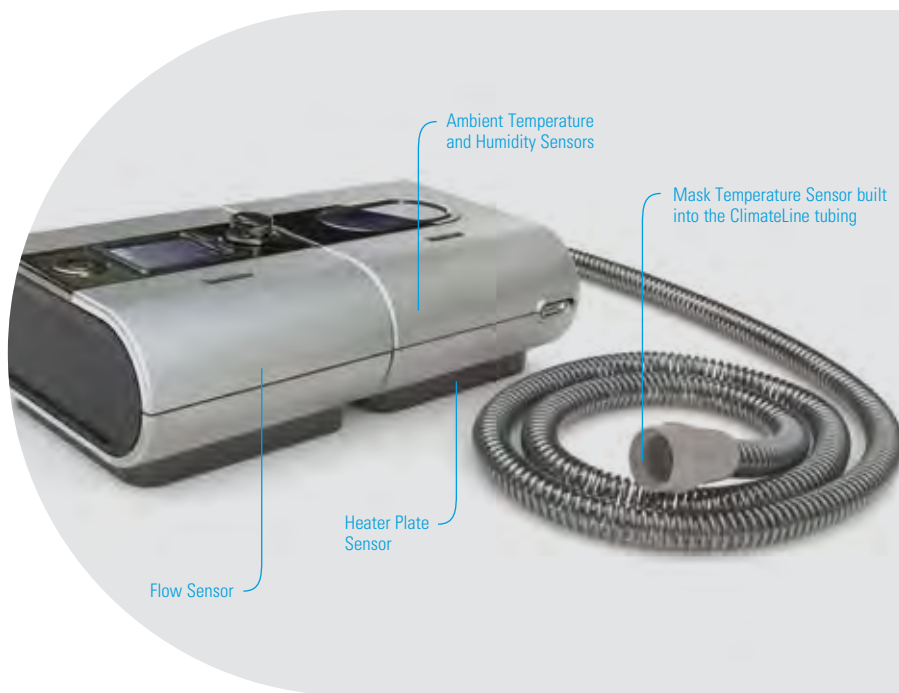

## Technology Highlight

### Climate Control **NEW**

Exclusive to ResMed, Climate Control prevents rainout without compromising humidity, providing the most comfortable therapy possible.

Technologies in the S9 therapy system and H5i™ heated humidifier work together with the new ClimateLine™ heated tube to automatically deliver the optimal temperature and humidity.

Unique temperature sensors ensure that the right conditions are maintained right to the mask, regardless of any ambient temperature and humidity changes.

## S9 Escape Auto—more value

The S9 Escape Auto is a whisper-quiet, elegant and comfortable therapy system for the treatment of Obstructive Sleep Apnea (OSA). The S9 Escape Auto provides Automatic Positive Airway Pressure (APAP) with Easy-Breathe Expiratory Pressure Relief (EPR) to prevent the collapse of your upper airway during sleep.

As the standard auto-adjusting model of the S9 Series, S9 Escape Auto offers excellent value with innovative features including Climate Control and the enhanced Easy-Breathe motor.

Together, these technologies deliver the S9 promise of more comfort.

## Specifications



<b>PERFORMANCE</b>	Operating pressure range: 4–20 cm H <sub>2</sub> O
<b>DATA STORAGE</b>	Compliance and therapy summary and statistic data: 365 sessions on device and SD Card
<b>DIMENSIONS (L X W X H)</b>	153 mm x 140 mm x 86 mm
<b>WEIGHT</b>	810 g
<b>90W POWER SUPPLY UNIT</b>	Input range 100–240V, 50–60Hz; 110V, 400Hz nominal for aircraft use; 70W (80VA) typical power consumption; 110W (120VA) maximum power consumption
<b>HOUSING CONSTRUCTION</b>	Flame retardant engineering thermoplastic
<b>OPERATING TEMPERATURE</b>	+5°C to +35°C
<b>AIR FILTER</b>	Standard: Polyester non-woven fiber Hypoallergenic: Acrylic and polypropylene fibers in a polypropylene carrier
<b>AIR TUBING</b>	SlimLine™ tube: Flexible plastic, 1.8 m, 15 mm inner diameter ClimateLine heated tube: Flexible plastic and electrical components, 2 m, 15 mm inner
<b>CLASSIFICATION</b>	Class II (double insulation), Type BF

## Product Codes: S9 Escape Auto

Australia, New Zealand	<b>36302</b>
Asia Pacific	<b>36362</b>